

SOFTCAIRE DX (DELUXE SYSTEM)

WEIGHT CAPACITY: 450 LBS



TABLE OF CONTENTS

IMPORTANT SAFEGUARDS	2
STATEMENTS AND SYMBOLS	3
INTRODUCTION/PACKAGE CONTENTS	4
PRODUCT FUNCTIONS	5
INSTALLATION	6
GROUNDING INSTRUCTIONS	7
OPERATION	8-10
• General	
• CPR Function	
• Pressure set-up	
• Low Pressure Warning	
CLEANING	11,12
HANDLING AND STORAGE	13
MAINTENANCE (General, Low Pressure)	13
TROUBLESHOOTING	14
WARRANTY	15

SAVE THESE INSTRUCTIONS

Please read and understand manual before using product.

IMPORTANT SAFEGUARDS

When using any electrical products, and especially when children are present, basic safety precautions should be followed. Those precautions should include the following:



DANGER

To reduce the risk of electrocution:

- Always unplug this product immediately after use.
- Do not use while bathing.
- Do not place or store this product where it can fall or be pulled into a tub or sink.
- Do not place or drop this product into water or other liquids.
- Do not reach for this product if it has fallen into water. Unplug immediately.

WARNING

To reduce the risk of burns, electrocution, fire or injury:

- This product should never be left unattended when plugged in.
- Close supervision is necessary when this product is used by, on or near children or invalids.
- Use this product only for its intended use as described in this manual.
- Do not use any attachments not recommended by the manufacturer.
- Never operate this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return the product to a service center for examination and repair.
- Keep the cord away from heated surfaces.
- Never block the air openings of the product or place it on a soft surface, such as a bed or couch, where the air openings may be blocked. Keep the air openings free of lint, hair and the like.
- Never drop or insert any object into any opening or hose.
- Do not use outdoors or operate where aerosol (spray) products are being used.
- Connect this product to a properly grounded outlet only. See Grounding Instructions.

GROUNDING

Before any connection is made, make certain that this product is connected to a protective earth conductor via the three-wire power cord; the three-blade plug should be inserted only into a socket outlet that provides a protective earth contact. Please refer to page 9 for more information.

STATEMENTS & SYMBOLS

NOTE, CAUTION, WARNING & DANGER STATEMENTS

NOTE

Indicates some tips or some information users should be aware of.

CAUTION

Indicate correct operating or maintenance procedure in order to prevent damage to or destruction of the equipment or other property.

WARNING/DANGER

Call attention to a potential danger that requires correct procedures or practices in order to prevent personal injury.

INTRODUCTION

This manual should be used for the initial set up of the system and for reference purposes.

INTENDED USE

SoftCaire DX pump and mattress is intended to reduce the incidence of pressure ulcers while optimizing patient comfort. It also provides the following:

- Individual home care setting and long-term care.
- Pain management as prescribed by a physician.

CONTRAINDICATION

Patient conditions for which the application of pressure relieving therapy on an alternation system is contraindicated are as follows:

- Cervical or skeletal traction
- Unstable spinal cord injuries

Manual: Please read and understand the manual before using this product.

PACKAGE CONTENTS

Check inside the package. If any of the following items should be missing or damaged, please contact your local dealer or service center for replacement immediately.

1. PUMP

SoftCaire DX pump and mattress is identical to which you purchased.



2. Mattress

The mattress comes with a zippered base. It has an 8" cell on cell construction with 5" upper layer alternating and 3" lower layer static.



3. Coversheet

Mattress cover is quilted and has a zipper so that it can be zippered to the base.



PRODUCT FUNCTIONS

PUMP

The functions of the SoftCaire DX pump are described below. Please refer to the below

figures.

(1) Power Switch

The switch is at the right side of the pump.

Turn ON/OFF the power and press Power button on the panel, the pump will start/stop operation.

(2) Press Cycle Time

The cycle time is selectable from the panel.

The cycle time value can be either 10, 15, 20, or 25 minutes.

(3) Press Pressure Range

20-55mmHg is offered from the panel, providing options of pressure value by easy press.

(4) Alarm Mute

The audible/visible alarm turns on when low pressure, power failure or alternate failure.

To mute the audible alarm, press the Alarm Mute button. The visible alarm indicator will flash till problem solved.

Re-press the Alarm Mute button to reactivate the alarm.

(5) Alternate mode

Press to set the air mattress in alternate mode.

(6) Static mode

Press to set the air mattress in static mode.

(7) Max Firm

Press to set the air mattress in quick inflation mode, which facilitates nursing and caring.

(8) Seat Inflate


Press to set the air mattress in seat inflate + alternate or seat inflate + static mode, pressure value will increase by 5mmHg, which facilitates seating.

Under max firm mode, keep max pressure of 55mmHg, no more 5mmHg increased.

(9) Lock button

Manual: press lock button for 3 seconds to lock the panel, press again for 3 seconds to unlock the panel.

1. Power switch
2. Pressure cycle time
3. Pressure range
4. Alarm mute
5. Alternate
6. Static
7. Max Firm
8. Seat inflation
9. Lock

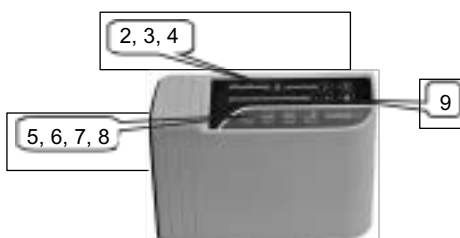

MATTRESS

SoftCaire DX pump includes a cell type mattress with materials made of (Nylon PVC, Nylon TPU or PU), provided optionally with the pump.

INSTALLATION

Step 1 Place the mattress flat on the bed frame. The inflation tube should be towards the foot end. The foot end is also indicated with a picture of footprints. See figure below.



Step 2 Cover the mattress system with a cotton sheet to avoid direct skin contact and reduce friction.

Step 3 Hang the pump unit on the foot board of the foot end of the bed frame and make sure it is stable.

Step 4 Connect the inflation tubes from the mattress to the pump's inflating nozzles. Make sure they are properly attached.

NOTE:

Make sure the air hoses are not kinked or tucked under the mattress. Also check to ensure the CPR valve is properly attached.

Step 5 Plug the power cord into an electrical outlet with grounded AC power. This product must be grounded. The power cord has a grounding wire with a grounding pin. This three-wire plug must be plugged properly into an outlet and grounded as shown in the figures in the Grounding Instructions section.

Important:

Before inserting the plug into the outlet, ensure the voltage is compatible.

Step 6 Turn on by pressing the power switch at the right side of the unit. Proceed to the Operation section.

Step 7 When not in use, disconnect the pump unit by unplugging the power cord.

GROUNDING INSTRUCTIONS

This product should be grounded. In the event of an electrical short circuit, grounding reduces the risk of electric shock by providing an escape wire for the electric current. This product is equipped with a cord that contains a grounding wire with a grounding plug.




DANGER Improper use of the grounding plug can result in a risk of electric shock.

If repair or replacement of the cord or plug is necessary, do not connect the grounding to either flat blade terminal. The wire with insulation having an outer surface that is green with or without yellow stripe is the grounding wire.

Check with a qualified electrician or serviceman if the grounding instructions are not completely understood, or if in doubt as to whether the product is properly grounded.

This product has a grounding plug that looks like the plug illustrated in figure A below. A temporary adapter, which looks like the adapter illustrated in figure B and C, may be used to connect this plug to a 2-pole receptacle as shown in figure B if a properly grounded outlet is not available. The temporary adapter should be used only until a properly grounded outlet (figure A) can be installed by a qualified electrician. The green colored rigid ear, lug, tab or the like extending from the adapter must be connected to a permanent ground such as properly grounded outlet box wire. Whenever the adapter is used, it must be held in place by the screw.

If it is necessary to use an extension cord, use only a 3-wire extension cord that has a three-blade grounding plug and a 3-slot receptacle that will accept the plug on the product. Replace or repair if damaged.

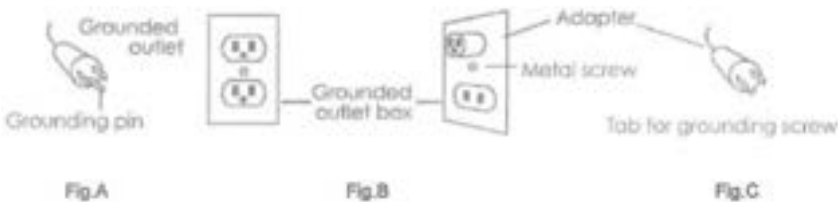

OPERATION (Always read the instructions before use.)

General

This product is designed to provide maximum comfort to home care patients. Please take care to operate this product properly to optimize its value.

Usage: Do not use another pump with different specifications. It is dangerous to use a pump with pressure capacity greater than 120mmHg. This may result in cell-puncture.

Do not change any component by yourself. If there is a need for replacement or repair, always contact your local dealer or service center.

For patients: When the Normal Pressure LED (green) comes on to indicate the pressure is adjusted to a desired firmness; the patient can then lie on the mattress.

Please cover the mattress with a cotton sheet to avoid direct skin contact and for the patient's comfort.

A firm surface will make it easier for the patient to transfer or reposition. Make use of the static mode function. It can support the patient by keeping the mattress from bottoming out when the patient is in a sitting position.

Hand check: Check to see if a suitable pressure is selected by sliding one hand between the air mattress and the foam base (or bed frame if there is no foam base) to feel the patient's buttock.

Users should be able to feel the space in between, and the acceptable range is approximately 25 to 40 mm (1" to 1-1/2"). This hand check procedure is issued by AHCPR.

OPERATION

Step 1

Turn on the power. A beep indicates it is beginning operation.

Step 2

The pump will enter static mode when the power is initially turned on. First it will inflate the mattress to 10-15mmHg and then switch to alternating mode.

If 10 mmHg cannot be reached within 30 minutes, the unit will stay in static mode and keep inflating. If 10-15 mmHg still cannot be reached within 40 minutes, the alarm will beep, and the LED will come on.

Press the Alarm Mute button to mute the alarm. The alarm LED will continue blinking. Press the Alarm Mute button again to reenable the audible alarm.

Step 3

The alternating cycle time and pressure level are preset to the middle settings. Select from the touch panel to adjust the cycle time and pressure level to the patient's specific requirements.

NOTE:

Press the static mode button from the touch panel to provide a firm surface that makes it easier for the patient to transfer or reposition. The static mode can support the patient by keeping the mattress from bottoming out when the patient is in a sitting position.

Step 4

During normal operation, the unit will monitor pressure. If the pressure is lower than the set pressure, the pump will automatically inflate the mattress to the set level and then stop. The alarm will beep, and its LED will come on to alert a low-pressure condition.

Press the Alarm Mute button to mute the alarm. The alarm LED will continue blinking. Press the Alarm Mute button again to re-enable the audible alarm.

Step 5

Press the Max Firm button to automatically inflate the mattress to the maximum level (above 55mmHg) for about 20 minutes. The pressure will return to a previously set level after 20 minutes.

In case there is a bad connection of tubing or any air leakage from the mattress, this will result in a low-pressure condition, and it will activate the alarm.

NOTE:

For suitable pressure, please refer to previous page for the hand check procedure.

OPERATION (CONTINUED)

CPR FUNCTION

When there is an emergency requirement to perform CPR on the patient, pull the two CPR straps at the head section of the mattress to release the air quickly from the mattress.

The CPR straps are located at the front right-hand side of the mattress.

PRESSURE SET UP

Users can adjust the pressure level of the air mattress to a desired firmness by themselves or according to the suggestion from a health care professional.

NOTE: It is recommended to press Auto Firm on the panel when the mattress is first inflated. Users can then easily adjust the air mattress to a desired firmness according to the patient's weight and comfort.

LOW PRESSURE WARNING

When abnormal pressure occurs, the Low-Pressure indicator will come on. For some types of pumps, the audible alarm will be activated to bring attention to a low-pressure condition.

Check if the connections are secure and correctly installed according to the relevant instructions.

NOTE: If the pressure is consistently low, open the zipper and confirm that all the hoses are properly connected. Then check for any noticeable leakage in any of the tubes. If necessary, contact your local dealer to replace any damaged tubes or hoses.

CLEANING

In this section, we describe the procedures to clean and decontaminate the pump. It is important to follow these procedures before using the system again. The cleaning task is required at least once a week to maintain personal hygiene.

PUMP

- DO NOT immerse or soak the pump in any water or fluids.
- Check for external damage and move the pump to the cleaning area.
- Place the pump on a work surface and spray or wipe the outside of the case with quaternary ammonium solution.
- DO NOT spray any cleaning solution directly on the surface of the pump.
- DO NOT use a Hypo carbonate or Phenolic based cleaning solution as this may cause damage to the case. Allow the solution to incubate for 10 minutes or accordingly as stated by the cleaning product instructions.
- Wipe case with a clean cloth. Make sure all areas are clean (top and bottom, both sides).
- Spray cloth with cleaning solution and clean faceplate. DO NOT allow excess cleaning solution on faceplate or control panel. (If solution gets inside, damage will occur.)
- Allow surface to thoroughly dry after cleaning.
- After the pump is thoroughly cleaned and dried, proceed to plug in the pump and test to see if it runs normally.
- Unplug the pump and store with proper identification tag.
- Avoid long exposure to sunlight.

MATTRESS

- Brush off or wipe down all surfaces of the cover sheet with soap and water before wetting with any liquid disinfectant.
- Any obvious blood spots should be wet thoroughly with 1:9 Hypochlorite solutions (1 part bleach to 9 parts water) and allow drying for at least 10 minutes. Then blot with a clean, damp cloth.
- Unzip the top cover from the mattress.
- Brush or wipe down all surfaces with soap and water before applying any liquid.
- Covers are immersed and soaked in disinfectant for the required incubation time.
- After pre-soaking, the cover is rinsed through a regular cycle in a washer with no soap then laundered with mild detergent (wash temperature 93°F, rinse temperature 78°F or on the coldest setting).
- Covers are aerated until they are fully dry. (Drying temperature range 90-120°F or on the coldest setting.)
- The air cells are unsnapped from one side and are sprayed on all sides with a disinfectant. Let it sit for the required incubation time and wipe down with a clean cloth. (Make sure to disconnect all the air cells, one by one, and spray the disinfectant on all sides, including all the connecting tubes and hoses. Let it sit for least 10 minutes.)
- If there is a base after you remove all the air cells, the base must be sprayed down with the disinfectant, inside and outside. Let it sit for the required incubation time and wipe down with a cloth.

Repeat the process with the tubing set: spray, incubate, and then wipe clean.

- The carrying bag should be turned inside out and completely wiped down using the disinfectant solution. Allow it to thoroughly air dry. Once the inside is dry, turn it back: wipe down the outside of the bag with disinfectant.

Dry the mattress on a flat surface area after cleaning, away from exposure to the sun. Avoid long exposure to sunlight)

HANDLING AND STORAGE

- Lay the mattress out flat and upside down.
- Roll from the foot end towards the head end; the foot-end strap can then be stretched around the rolled mattress to prevent unrolling.
- Do not fold, crease or stack the mattress.

MAINTENANCE

General

- Check the power cord and plug to see if there are abrasions or excessive wear.
- Check the mattress cover for signs of wear or damage. Ensure the mattress cover and tubes are connected correctly.
- Plug in the pump unit and check the airflow from the hose connection port. The airflow should alternate between ports every half-cycle time.
- Check the air hoses to see if there are any kinks or breaks. For replacement, please contact your local agent or dealer.
- Make sure the mattress tube is well connected.
- Check the pump unit and make sure both power indicators are off when the switch is turned off.

Low pressure

Examine if there is any air leakage between the pump and the mattress connections or from the air mattress tubes:

- Check connectors between the air mattress and pump. If there is any disconnection, please reconnect it.
- Check the CPR Valves. Ensure their outlets are sealed.
- Check the air-connecting tubes. Ensure each single cell is not broken.
- Set the pressure-selector knob to Firm position. Keep the tubes fully inflated and inspect for air leakage.

Check for air leakage from cells. Ensure no leakage occurs. If any leakage occurs, please contact your local agent or dealer.

TROUBLESHOOTING

Problem	Solution
The pump doesn't work	<ol style="list-style-type: none">1. Ensure the plug is inserted firmly into the outlet.2. Turn on the power switch again. <ul style="list-style-type: none">• If the power LED is ON and the pump still doesn't work, contact your local dealer immediately.• If the power LED is OFF, there may be a faulty outlet. Try to connect the power cord to another outlet. If the power LED is still OFF, contact a qualified electrician for main power check.
Incomplete Inflation (Low Pressure)	<ol style="list-style-type: none">1. For a quick check, adjust the pressure to Firm position.2. Check to see if the tubes connected to the pump are twisted or there is any leakage occurring. <ul style="list-style-type: none">• Always keep the tubes straight.• Change tubes if there is any leakage.• Ensure the CPR valves are closed. Ensure every single cell is not broken.
Slow air flow	<p>A dirty filter may decrease the air flow. Wash the filter with mild detergent to keep it clean.</p> <ul style="list-style-type: none">• Check the air filter at the back of the pump unit at least once a month.

WARRANTY

Your Rhythm product is warranted to be free of defects in materials and workmanship for eighteen (18) months to the original consumer purchaser. This warranty does not cover device failure due to misuse, negligence, improper storage or handling, improper operation, and unauthorized modifications. This warranty does not cover device failure due to owner misuse or negligence, or normal wear and tear.

The warranty does not extend to non-durable components, such as rubber accessories, which are subject to normal wear and need periodic replacement.

For warranty service, it is recommended but not required, that the product be returned to the Rhythm Healthcare dealer through whom it was originally purchased. However, all Rhythm Healthcare dealers are qualified to assist you in obtaining warranty service. If the product is to be returned to Rhythm Healthcare, prior authorization will be required.

This product is identified by a serial number at the time of manufacture. Removing or defacing the serial number label may void this warranty. If you have a question about your Rhythm Healthcare product or this warranty, please email us at contactus@rhythmhc.com



Rhythm Healthcare
Clearwater, FL 33762

Email: contactus@rhythmhc.com
Website: www.rhythmhc.com

V2: 7/25/23